



Gooch & Housego



450 – 750nm Acousto-Optic Tunable Filter

TF600-300-6-3-GH19A

An acousto-optic tunable filter for use at 450 – 750nm. Ideal for wavelength selection of Supercontinuum sources, particularly in confocal microscopy applications.

Designed to optimise diffracted beam efficiency, quality & pointing stability while maintaining high extinction ratio and low RF drive power, also incorporating patented side-lobe suppression and active temperature stabilisation.

Use in conjunction with the 64040-150-10MDFS-16X1 frequency synthesised driver, enables independent control of frequency, amplitude and phase of up to sixteen channels, allowing active pass band resolution and profile control.

Gooch & Housego's AOTF capability is extensive. By combining our scientific knowledge, modelling capability and engineering expertise with our renowned manufacturing skill and high quality, our products are aimed at the most discerning customers, in the most demanding applications.

Key Features:

- 450 – 750nm
- High Efficiency
- Excellent pointing stability
- Patented side-lobe suppression
- High extinction ratio
- Maintains beam quality
- Low drive power

Applications:

Biomedical:

- Confocal Microscopy

General Specifications


Interaction material:	Tellurium Dioxide (Anisotropic)
Wavelength range:	450 - 750nm
Resolution (FWHM):	< 6nm at 525nm
Active aperture:	3mm
Incident polarisation:	Linear, vertical with respect to base
Polarisation of diffracted order:	Linear, orthogonal to input (90° rotated)
Pointing stability of diffracted order:	< ±0.01°
Beam separation:	> 3°
Diffraction efficiency:	> 90%
RF drive power:	< 100mW / channel

Recommended RF Driver: 64040-150-10MDFS-16X1

Ordering Code

Explanation: TF600-300-6-3-GH19 (AOTF, wavelength range 450 – 750nm, 6nm resolution, 3.0mm active aperture, GH19 housing, side lobe suppression).

T	F	6	0	0	-	3	0	0	-	6	-	3	-	G	H	1	9	A
---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---	---

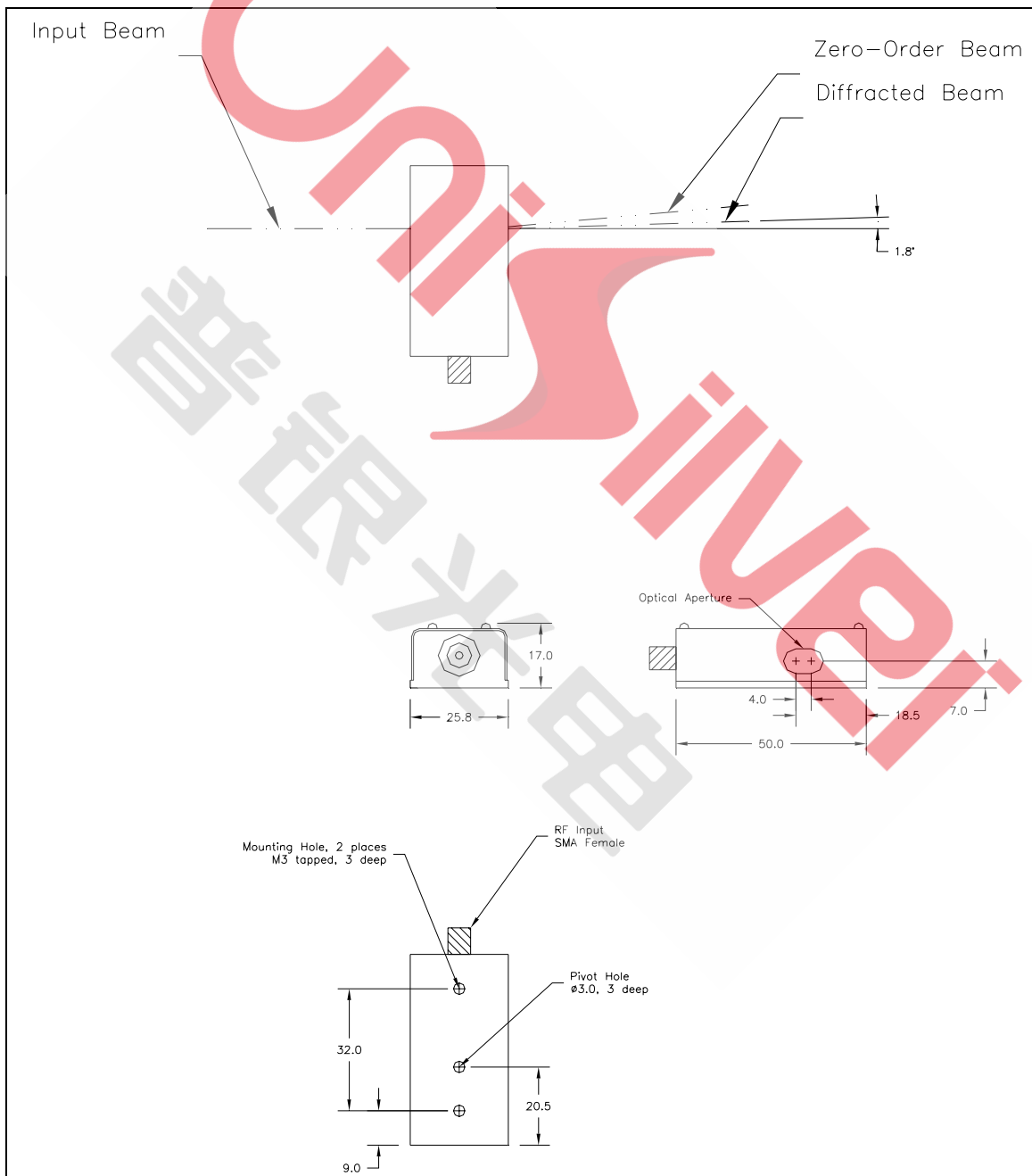

Tuning relation

